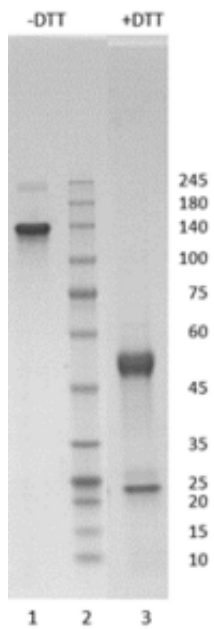




## Rabbit IgG1-kappa antibody to SARS CoV-2 NP (clone 82C3)

Catalogue #	R1-166-100
Immunogen:	Nucleocapsid protein SARS-CoV-2
Immunogen Description:	Recombinant SARS-CoV-2 Nucleocapsid protein produced by CHO-based Icosagen Cell factory Ltd. proprietary suspension cell line
Clonality:	Rabbit monoclonal
Clone:	82C3
Class:	rlgG1
Reactivity:	SARS-CoV-2 Nucleocapsid protein
Dissociation constant (K <sub>D</sub> ):	<1.0 x 10 <sup>-12</sup> M
Application:	ELISA
Kd:	<1.0 x 10 <sup>-12</sup> M
ELISA:	indirect ELISA 0,02-0,04 ng/ml
Purification:	Protein A affinity chromatography following desalting
Buffer:	PBS pH 7.4
Shipping:	Shipped at ambient temperature.
Storage:	Store at +4 °C. Avoid multiple freeze-thaw cycles.

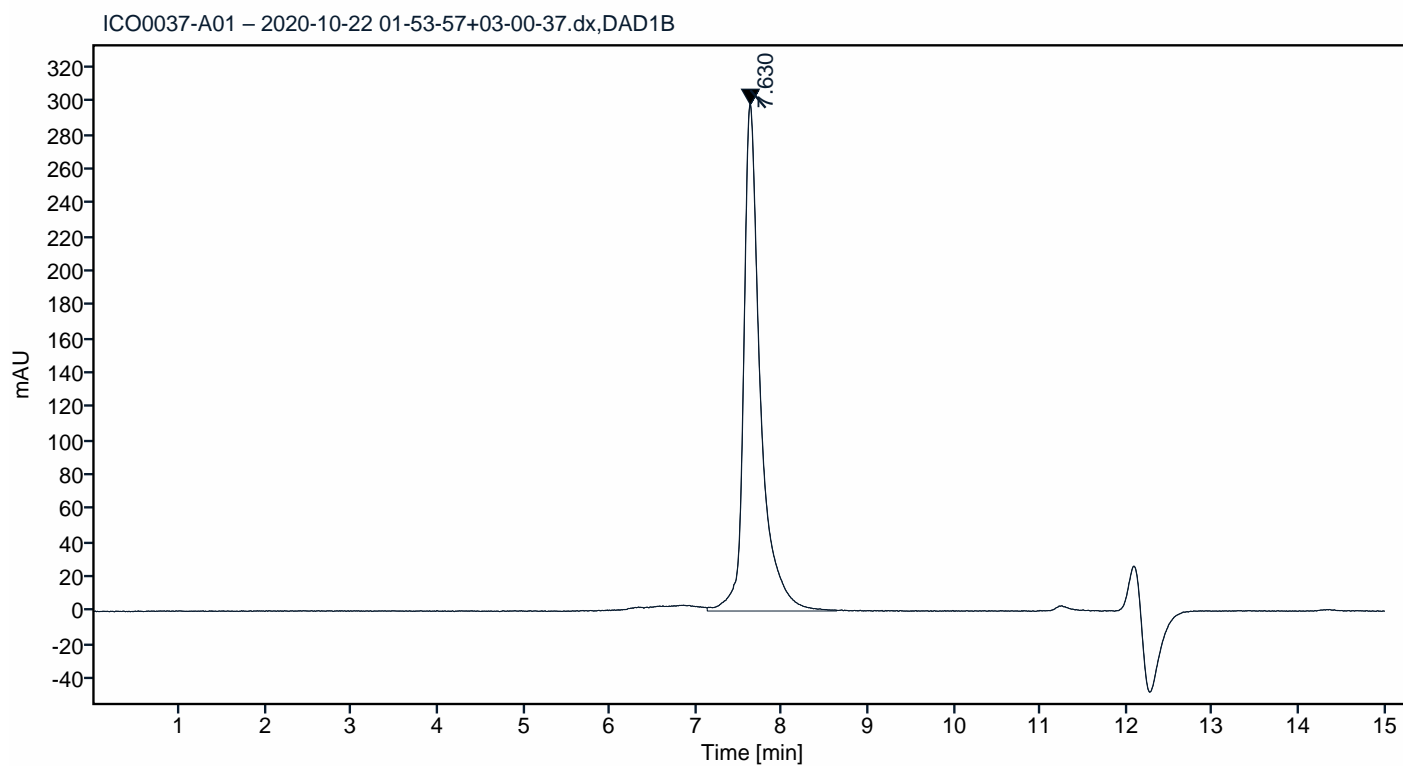


**Figure 1.** Simply Blue Safe stained SDS-PAGE analysis of monoclonal antibody to SARS-CoV-2 Nucleocapsid protein, clone 82C3. 4-12% gradient gel is used for analysis. Lane 1. Monoclonal antibody to SARS-CoV-2 Nucleocapsid protein, clone 82C3 (-DTT). Lane 2. Size marker. Lane 3. Monoclonal antibody to SARS-CoV-2 Nucleocapsid protein, clone 82C3 (+DTT).

## Peak Table

Peak #	RT (min)	Area	Area %
1	7.630	4033.86	100.00

## Chromatogram

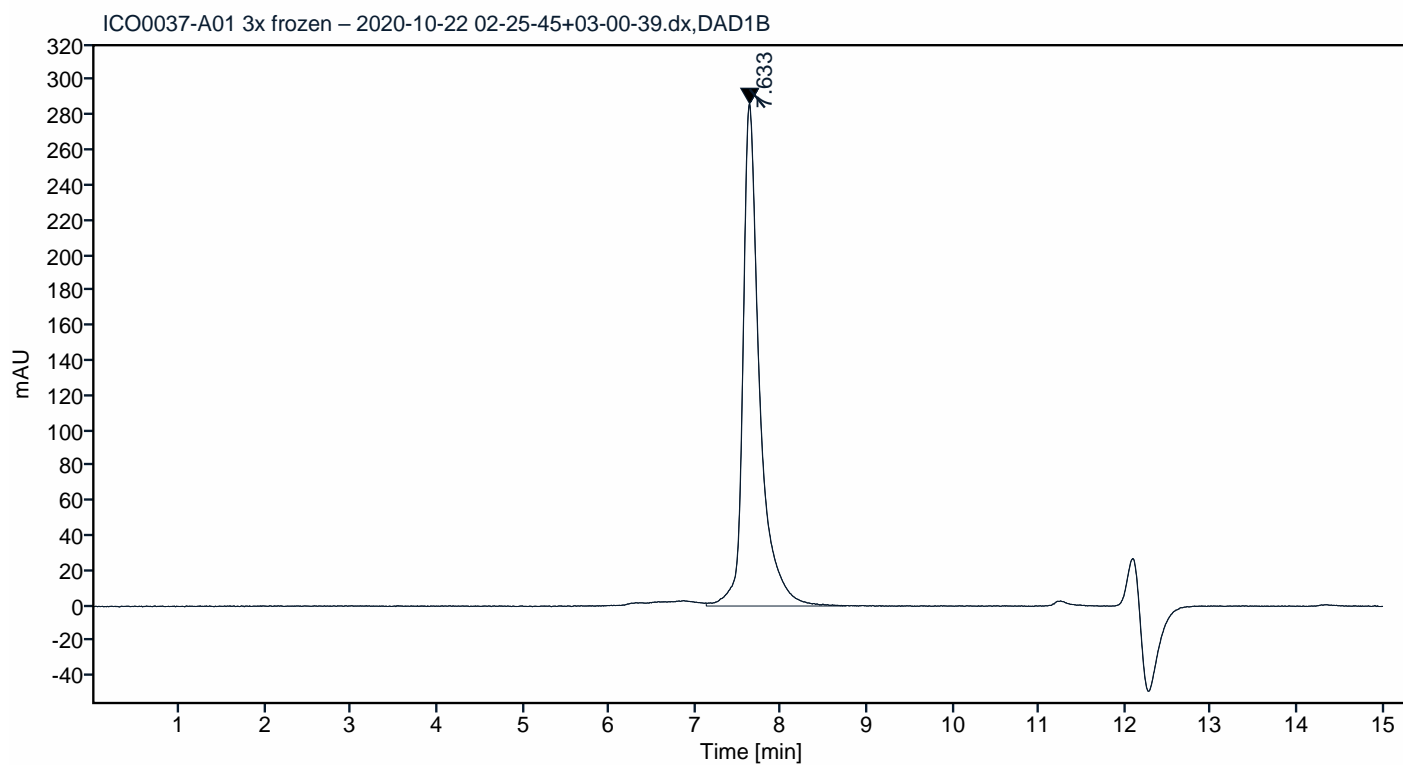


**Figure 2.** HPLC analytical SEC for final product.

## Peak Table

Peak #	RT (min)	Area	Area %
1	7.633	3867.96	100.00

## Chromatogram



**Figure 3.** HPLC analytical SEC after 3 freeze-thaw cycles.

**Figure 4.** Octet RED96e analysis. Antibody was loaded on sensor for capture of Nucleocapsid protein.